
ID DL 1766.LT.ALB.EAD`S.2015

PACHYMETER SP-100



- Automatic CCT corrected IOP calculations
- Wide measuring range 150 to 1,200 µm
- Patient data input & export capability
- Built-in printer
- 8 readings with averaging
- Cordless

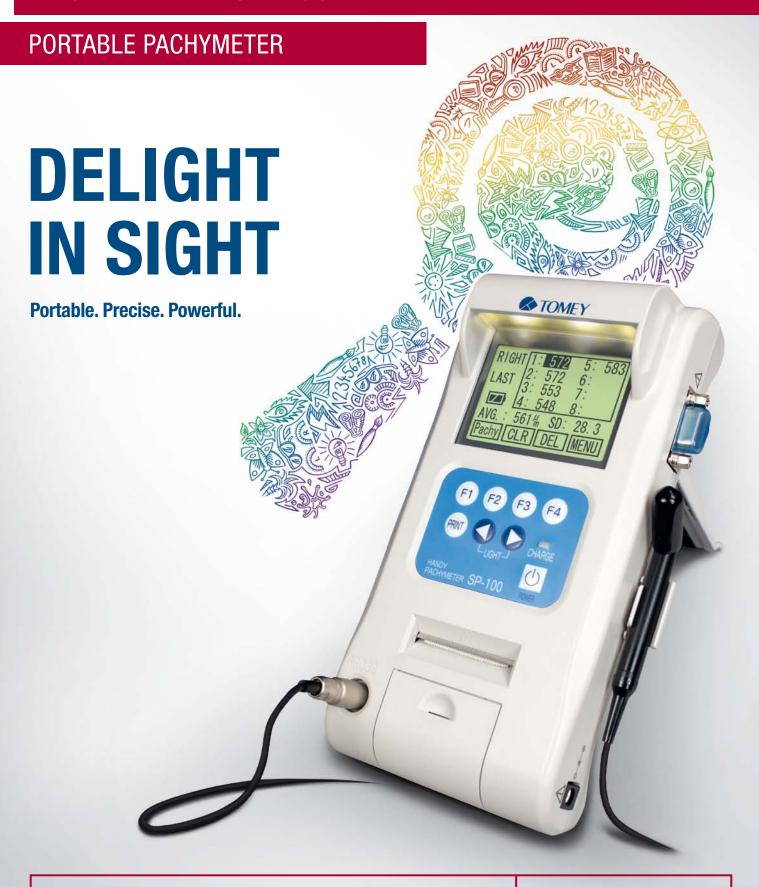

THE TOMEY PORTABLE PACHYMETER SP-100



QUALITY IN DETAIL

Small Size - High Ability


The Tomey pachymeter **SP-100** covers all your needs regarding pachymetry. Small and portable design, wireless operation and light weight make the **SP-100** an "easy-to-use unit".

Customised Pachymeter


You need adjustable IOP correction? Easily done! With the programmable IOP formula you are able to get the central corneal thickness (CCT) corrected IOP. The **SP-100** calculates this automatically with the formula of your choice. The system displays the average of 8 measurements. Every single measurement is already an average out of up to 20 exams.

Documentation


You can print your results on the built-in printer or transmit it to a PC or network system via serial port. The printout includes the patient ID, name and measured data. The rechargeable battery with AC/DC adapter has a standby time up to 6 hours.

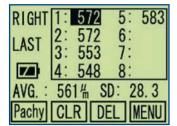

Average will be calculated out of 8 single measurements

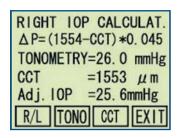

IOP calculation

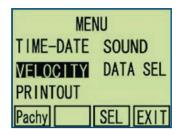

Settings display



Illuminated LC display

2004/ 9/24 16:30 ID:12345678901234 NAME **ABCDEFGHIJKLMN** Vel:1640 Average RIGHT 1: 548 2: 572 3: 553 4: 548 5: 583 Avg: 561 SD:28.3 LEFT -1: 823 2: 830 3: 826 Avg: 826 SD: 2.9

Sample printout



Control panel



Pachymeter probe

SPECIFICATIONS

PACHYMETRY FUNCTION

Single Mode 8 single, average / minimum / latest

minimum / latest options 20 MHz

 Frequency
 20 MHz

 Velocity range
 1,400 to 2,000 m/s

 Default setting
 1.640 m/s

Measurement range1,400 to 2,000 m/sMeasurement accuracy± 0.005 mm

MAIN UNIT

Display STN 2.7" dot Matrix LCD 240 x 160 pixel

PrinterBuilt-in thermal printerCommunication portRS-232

DIMENSIONS & ELECTRIC REQUIREMENTS

Dimensions WDH118 x 207 x 59 mmWeightapprox. 0.53 kgBatteryModel Sonyo GSWorking time6 hours 45 minutes

Battery charger

Output7.0 VLoad current2A

TOMEY EUROPE TOMEY GmbH Am Weichselgarten 19a 91058 Erlangen, Germany

120

TOMEY ASIA-PACIFIC TOMEY CORPORATION JAPAN 2-11-33 Noritakeshinmachi Nishi-ku, Nagoya 451-0051, Japan

Phone +81 52 581 5327, Fax +81 52 561 4735 Email intl@tomey.co.jp

